

Mount-Quick Aqueous

Catalog No: tcdt312



Available Sizes

Size: 1Each



Specifications

Application:

IHC Slide Mounting

Form:

Viscous Gel

Storage Instruction:

After use, always replace cap and store in a cool, dark place.

Protocol:

Drop onto the cover glass or section from above. (Prevent bubble formation by inverting the container and touching the tip onto the bottom of the drop.) * Before first use, cut open the container tip at a 45° angle with a box cutter or scissors.

Notes

Mount-Quick Aqueous can be used with methyl green and BCIP/NBT staining, but are not suited for long-term storage. After staining, we recommend photographing as soon as possible, as the color can fade in 24 hours after drying, dehydrating, clearing, and mounting. BCIP/NBT can stain red and black with immunostaining.

Product Description


Mount-Quick Aqueous is Sample mounter for permanent preservation.

(water soluble mounting medium for mounting after staining without dehydration)

Maintains superior transparency when fully hardened. Low fluorescence.

Mounting is simple with easy-to-use dropper bottle.

Specifications: Mobility: pale yellow, transparent, viscous aqueous solution; specific gravity: 0.88sp; viscosity: 33cp; refractive index: 1.41

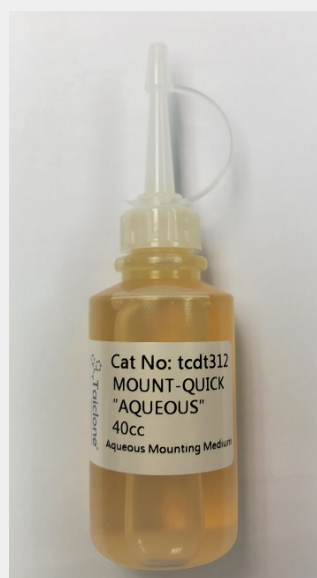

All products are for RESEARCH USE ONLY. Not for diagnostic & therapeutic purposes!